


Prelab 2

- | Reagent Name | Reqd. PPE | Critical Safety Hazards | Reactivity | Disposal |
|--------------|--|-------------------------|-------------|--|
| FD&C Red 40 | Standard PPE
(gloves, eyewear,
lab coat) | None | Nonreactive | Can be disposed of
in sink since it is
food-safe |

2. Dilution Procedure

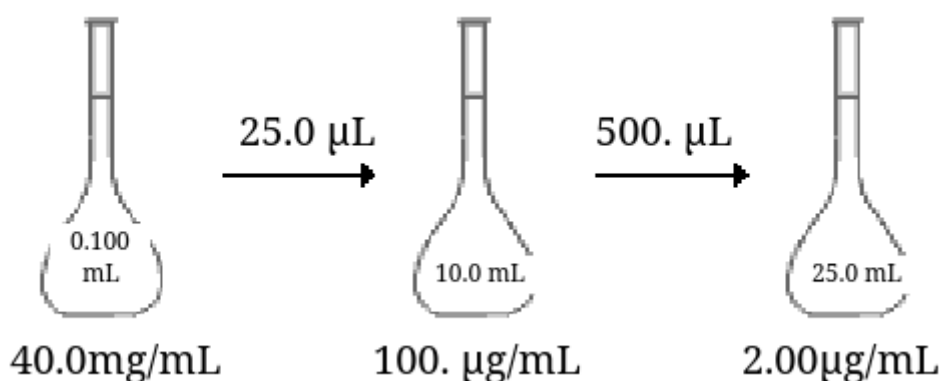

$$40.0 \text{ mg/mL} \times \frac{25.0 \mu\text{L}}{10.0 \text{ mL}} \times \frac{500. \mu\text{L}}{25.0 \text{ mL}} = 2.00 \times 10^3 \text{ mg/mL} \times \left(\frac{\mu\text{L}}{\text{mL}} \right)^2 \times \left(\frac{1 \text{ mL}}{1000 \mu\text{L}} \right)^2 \times \frac{10^6 \text{ ng}}{1 \text{ mg}} = 2.00 \times 10^3 \text{ ng/mL}$$

- Fill a 50 mL beaker with around 40 mL of DI water; set aside.
 - Using a P50 micropipette, dispense 25.0 μL of the 40.0 mg/mL food coloring solution into a 10.0 mL volumetric flask.
 - Using the DI water from Step 1, fill that volumetric flask to the mark. A transfer pipette should be used when approaching the mark for finer control to avoid overshooting.
 - Cap the volumetric flask and invert 20 times to mix.
 - Using the P1000 micropipette, transfer 500 μL from the volumetric flask into a new 25.0 mL volumetric flask. Using the DI water from Step 1, fill that volumetric flask to the mark, using a transfer pipette when approaching the mark to avoid overshooting.
 - Cap the volumetric flask and invert 20 times to mix.
 - Four cuvettes should be prepared for UV-Vis by wiping the clear sides of the cuvette with a KimWipe and ensuring that only the frosted sides are touched from that point until retrieval from the UV-Vis machines.
 - Fill a cuvette with DI water and use it as the blank to calibrate the UV-vis
 - Use a transfer pipette to fill each remaining cuvette with about 3 mL (3 cm) of solution.
 - Use the UV-Vis machine to record the absorption for each cuvette at 503 nm.
3. $A_\lambda = c\varepsilon_\lambda l$
- $$A_{503 \text{ nm}} = (2.00 \times 10^3 \text{ ng/mL})(2.59 \times 10^4 \text{ L mol}^{-1} \text{ cm}^{-1})(1.0 \text{ cm})$$
- $$A_{503 \text{ nm}} = (2.00 \times 10^3 \text{ ng/mL})(10^{-9} \text{ g/ng})(10^3 \text{ mL/L})(2.59 \times 10^4 \text{ L mol}^{-1} \text{ cm}^{-1})(1.0 \text{ cm})$$
- $$A_{503 \text{ nm}} = (2.00 \times 10^{-3} \text{ g/L})(2.59 \times 10^4 \text{ L mol}^{-1} \text{ cm}^{-1})(1.0 \text{ cm})$$
- $$A_{503 \text{ nm}} = (20.0 \text{ g})(2.59 \text{ mol}^{-1})(1.0)$$
- $$A_{503 \text{ nm}} = (20.0 \text{ g})(2.59 \text{ mol}^{-1})(1.0) \frac{1 \text{ mol}}{496.42 \text{ g}}$$
- $$A_{503 \text{ nm}} = 0.10$$
4. The final solution of my SOP has a concentration of $2 \times 10^3 \text{ ng/mL}$. This is the same as in Problem 3. Therefore, I expect the solution to have an absorbance of 0.10.

5. My answers are the same because the questions are the same.